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The RETSCH AS 400 control is used for sieving dry goods with test sieves up to 400 mm in diameter. The uniform, horizontal circular sieving motion produces a sharp separation of the sample fractions. Fine and coarse-grained goods from areas such as milling, brewing, chemical industry, quarries, soil testing, woodworking and plastics industry, can be exactly separated with the AS 400 control. This particular sieving motion is preferably used for long or fibrous, needle-shaped or flat materials due to their horizontal orientation. For the testing of plastics (grainy molding materials), the standard DIN 53 477 stipulates exactly this circular sieving motion.

The AS 400 control can be used as test instrument for the quality control according to DIN EN ISO 9001. Due to the controlled drive which is independent of the power frequency the AS 400 control yields reproducible results worldwide. The speed and sieving time are set, displayed and monitored digitally. The instrument is supplied with a test certificate and can be recalibrated. If desired, the rotation direction can be set to alternate in the interval. A memory for 9 sieving programs facilitates routine analyses. The AS 400 control has an integrated interface for controlling all sieving parameters via the EasySieve® software.

ACCURACY & EFFICIENCY

- | Circular sieving motion according to DIN 53477
- | For sieves up to (Ø) 400 mm
- | Measuring range 45 µm to 63 mm
- | Easy operation, ergonomic design
- | Free digital selection of process parameters (time, speed, interval)
- | Memory for 9 Standard Operating Procedures (SOPs)
- | Test materials monitoring according to DIN EN ISO 9001

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IQ/OQ Documents

We provide IQ/OQ documentation for the „control“ sieve shakers to support IQ/OQ certification by our customers.

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Ultrasonic baths and dryers


Suitable for thorough cleaning of test sieves and for quick, gentle drying of samples and sieves.

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RETSCH TEST SIEVES AND ACCESSORIES - ENGINEERED FOR SUPERIOR PERFORMANCE

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
EasySieve, the software for particle size analyses, exceeds manual evaluation in many aspects. The software is able to automatically control the necessary measurement and weighing procedures – from the registration of the weight of the sieve up to the evaluation of the data. It is simple and convenient to use and is also available in an FDA 21 CFR Part 11-conform version.

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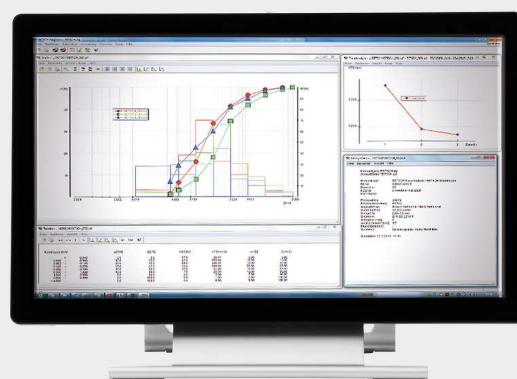
RETSCH's Analytical Sieve Shaker AS 400 control is ideally suited for separation, fractioning, particle size determination of construction materials, ground corn, plastics, sand, splints

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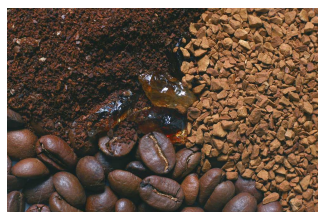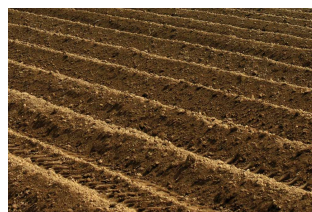To find the best solution for your sample preparation task, visit our application database.

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HORIZONTAL SIEVE SHAKER AS 400 CONTROL

Horizontal Sieve Shaker AS 400 control for test sieves up to 400 mm Ø
(please order clamping device, test sieves and collecting pan separately)

30.022.0001

AS 400 control

100–240 V, 50/60 Hz

incl. test report acc. to EN
10204 2.2

CLAMPING DEVICES AS 400 CONTROL

32.662.0010

Clamping device "standard", for test sieves 400 mm Ø

32.662.0011

Clamping device "comfort", for test sieves 400 mm Ø

SIEVE STACKS AND ACCESSORIES AS 400 CONTROL

60.166.000999

Sieve stack consisting of 6 test sieves (ISO 3310-1), 400 mm Ø, 65 mm height
(500 µm, 1 mm, 2 mm, 5 mm, 10 mm, 20 mm) and collecting pan

60.167.000999

Sieve stack consisting of 6 test sieves (ASTM E11), 400 mm Ø, 65 mm height
(35 mesh, 18 mesh, 10 mesh, 4 mesh, 3/8", 3/4") and collecting pan

ACCESSORIES AS 400 CONTROL

03.802.0006

Transportation handles AS 400

99.200.0026

IQ/OQ Documentation for AS 400 control

ACCESSORIES FOR CLAMPING DEVICES AS 200, AS 300, AS 400

CLAMPING LIDS

32.481.0022

Clamping lid with large window of Perspex for test sieves 200/203 mm Ø

02.660.0012

Clamping lid with large window of Perspex for test sieves 400 mm Ø

CLAMPING ELEMENTS

32.142.0001		Clamping nuts, (2 pieces) for clamping device "standard"
32.737.0001		Quick-clamping elements, (2 pieces) for clamping device "comfort" AS 200/300/400
05.114.0080		O-ring for quick-clamping element for AS 300 and AS 400, 1 piece

CLAMPING RODS

32.248.0002		Threaded rods, (2 pieces) for clamping device "standard"
32.248.0001		Threaded rods, short, (2 pieces) for clamping of max. 5 test sieves for clamping device "standard"
32.742.0010		Rods, smooth, (2 pieces) for clamping device "comfort" AS 300/400

SIEVING AIDS

32.365.0001		Chain ring for test sieves 200 mm and 203 mm Ø to support horizontal sieving
32.050.0001		Brushes, 3 pieces
32.902.0001		Cubes of polyurethane, 12 x 12 x 12 mm, 10 pieces
32.902.0002		Cubes of polyurethane, 20 x 20 x 20 mm, 10 pieces
32.354.0001		Balls of rubber, 20 mm Ø, 5 pieces
32.354.0002		Balls of agate, 10 mm Ø, 10 pieces
32.354.0004		Balls of steatite, 6 mm Ø, 150 g

ACCESSORIES FOR TEST SIEVES (PANS, RINGS, LIDS)

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