




КРИОГЕННА МЕЛНИЦА

Efficient cryogenic grinding at -196°C

Cryogenic grinding is a process where thermally sensitive and elastic substances are successfully processed by cooling with liquid nitrogen. The CryoMill is a laboratory ball mill specifically designed for this application. It features an integrated cooling system which continually cools the grinding jar with liquid nitrogen before and during the grinding process. Thus, the sample is embrittled and volatile components are preserved.


The liquid nitrogen is continually supplied from an autofill system in the exact amount required to keep the temperature at $-196\text{ }^{\circ}\text{C}$. The user never comes into direct contact with LN_2 which ensures a high degree of operational safety. The CryoMill's versatility (cryogenic, but also wet and dry grinding at room temperature) makes it the ideal grinder for sample quantities up to 20 ml. Powerful impact ball milling results in highest grinding efficiency.

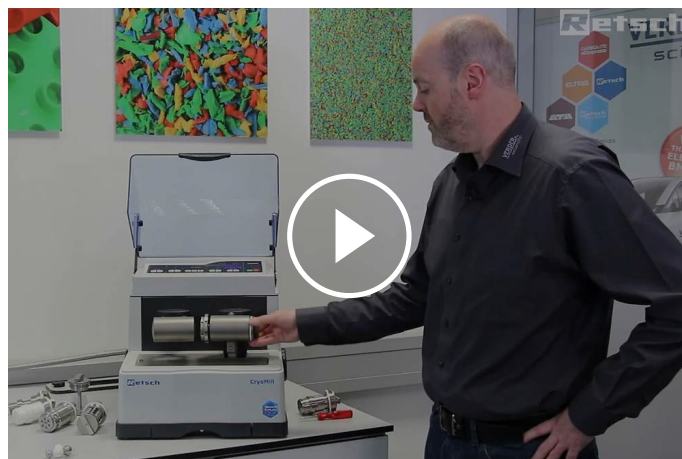

[Click to view video](#)

UNMATCHED PERFORMANCE FOR CRYOGENIC BALL MILLING

- | Max. speed 30 Hz
- | Horizontal oscillation causes strong impact effects for effective sample processing
- | Up to 8 mm feed size and $5\text{ }\mu\text{m}$ final fineness
- | 1 grinding station for jars of min. 5 ml and max. 50 ml, adapter for $6 \times 2\text{ ml}$ tubes
- | Jars are continuously cooled with liquid nitrogen during processing, safe handling thanks to autofill tank
- | Small bench top model
- | Storable SOPs and cycle programs, compact bench top model, 4 different jar materials for dry and wet grinding

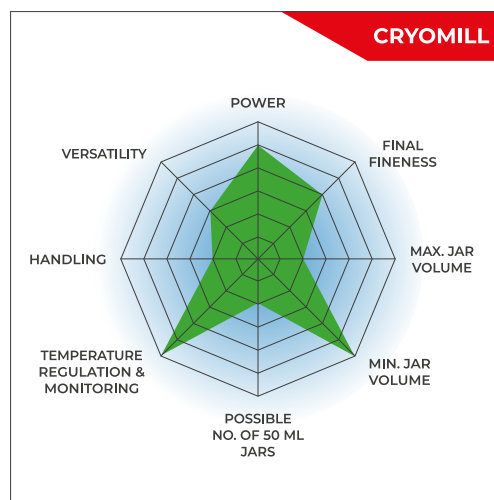

VERSATILE & SAFE

- | powerful cryogenic grinding by impact and friction with up to 30 Hz
- | 3 different grinding modes (cryogenic, dry/ wet at ambient temperature)
- | closed LN₂-system (Autofill) for enhanced user safety
- | screw-top jars ensure convenient, leak-proof operation
- | wide range of accessories including various jar and ball sizes, adapter racks and LN₂ feeding system
- | optional zirconium oxide jars specifically designed for cryogenic grinding
- | 9 programmable cooling and grinding cycles (10 s to 99 min)



USER SAFETY & PERFECT RESULTS GUARANTEED

The CryoMill is a cryogenic grinder designed with user safety in mind. The liquid nitrogen flows through the closed system and the user never comes into direct contact with LN₂ which ensures a high degree of operational safety. The automatic cooling system guarantees that the grinding process is not started before the sample is thoroughly cooled. This results in reduced consumption and guarantees reproducible cryogenic grinding results.

The CryoMill is very easy to operate. Parameters such as oscillation frequency, pre-cooling or grinding time can be digitally set via a clearly structured keypad. LEDs in the display indicate the current state of operation, e.g. cooling or grinding. Usually, grinding only takes a few minutes so that the sample does not get warm during the process. If, however, longer grinding times are required, it is also possible to pre-select periods of intermediate cooling and the number of cryogenic cycles.



[Click to view video](#)

All instrument parameters are retained during standby operation for subsequent processes. The laboratory mill can also be operated without cooling which makes it suitable for a vast range of applications.

FOR SAFE AND EFFECTIVE GRINDING PROCESSES
ACCESSORIES FOR THE CRYOMILL



GRINDING JARS & ADAPTERS

The CryoMill is equipped with one grinding station for screw-top grinding jars with volumes of 10 ml, 25 ml, 35 ml or 50 ml. It is also possible to use adapters for 4 grinding jars of 5 ml each as well as for 6 reaction vials of 2 ml each. A 25 ml grinding jar of zirconium oxide and matching grinding balls as well as a PTFE jar are available for applications where steel would cause sample contamination.



LIQUID NITROGEN FEED

For safe and comfortable operation, RETSCH provides an autofill system for liquid nitrogen which is available with a 50 liter container and provides cooling in cryogenic grinding applications for approximately 5 hours. It is also possible to connect existing cryo tanks to the mill, using a connection tube with safety valve.

КРИОГЕННА МЕЛНИЦА

ТИПИЧНИ МАТЕРИАЛИ ЗА ПРОБА


Due to the automatic embrittlement of the samples during cryogenic grinding the CryoMill is suitable for pulverizing, for example, waste, soil, chemical products, tissue, hair, wood, sewage sludge, bones, plastics, oil seed, paper, plants, pills, textiles, animal feed, wool etc.



rubber duck



hard plastic

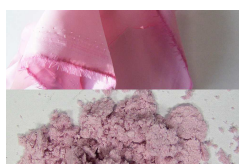

*текстилни
изделия*



gummy bears



Parsley

[Click to view video](#)


To find the best solution for your sample preparation task, visit our application database.

КРИОГЕННА МЕЛНИЦА

ПРИНЦИП НА РАБОТА

The grinding jar of the CryoMill performs radial oscillations in a horizontal position. The inertia of the grinding balls causes them to impact with high energy on the sample material at the rounded ends of the jar and pulverize it.

The grinding jar is continually cooled with liquid nitrogen from the integrated cooling system before and during the cryogenic grinding process.

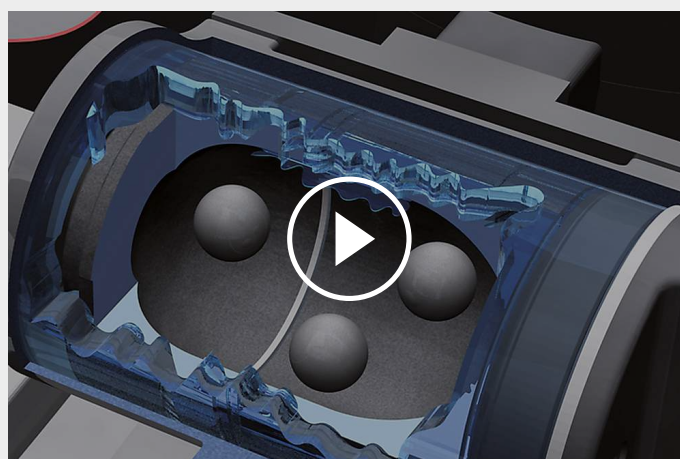

[Click to view video](#)

КРИОГЕННА МЕЛНИЦА

TECHNICAL DATA

Приложения	намаляване на размера, смесване, хомогенизиране, разрушаване на клетки
Област на приложение	биология, геология / металургия, инженеринг / електроника, медицина / фармация, околната среда / рециклиране, селското стопанство, строителни материали, стъкло / керамика, химия / пластмаси, храна
Материали	твърд, средно твърди, мек, крехък, еластични, влакнести
Принцип на смилане	удар, триене
Максимален размер на захранващият продукт	<= 8 мм
Големина след смилане*	~ 5 микрона
Размер на партидата / количество на пробата*	макс. 20 мл
Брой на смилачните станции	1
Vibrational frequency	дигитално, 5 - 30 Hz (300 - 1800 об/мин)
Най-често използвано време за смилане	10 мин / 4 мин (охлаждане / смилане)
Сухо смилане	да
Мокро смилане	да
Криогенно смилане	да
Разрушаване на клетки	да
Защита на затягащото устройство	да
Видове смилачни цилиндри	с резба
Материал на смилачните части	закалена стомана, неръждаема стомана, циркониев оксид, PTFE
Обем на смилачната камера	5 мл / 10 мл / 25 мл / 35 мл / 50 мл
Автоматична система за пълнене	50 л
Настройване времето за смилане	дигитално, 30 сек - 99 мин
Запаметяване на настройки	9
Данни за електрозахранването	100-240 V, 50/60 Hz
Захранване	еднофазово
Защита	IP 30
Енергийна консумация	260 W

Размери затворена	395 x 373 x 577 мм (D: 710 мм с изпускателна тръба)
Нетно тегло	~ 45 кг
Стандарти	CE

в зависимост от захранващият материал и настройките на инструмента

www.retsch.com/cryomill

ДАННИ ПОРЪЧКА

CRYOMILL

(please order Autofill with LN2 container and safety valve, grinding jars and balls separately)

20.749.0001



CryoMill, 100–240 V, 50/60 Hz

GRINDING JARS CRYOMILL

HARDENED STEEL

01.462.0300



5 ml, to be used with adapter 02.706.0304

01.462.0330



25 ml

01.462.0329



35 ml

01.462.0328



50 ml

STAINLESS STEEL

01.462.0290

5 ml, to be used with adapter 02.706.0304

01.462.0331



10 ml

01.462.0334



25 ml

01.462.0333



35 ml

01.462.0332



50 ml

ZIRCONIUM OXIDE

01.462.0336



25 ml

PTFE

01.462.0335



25 ml

ACCESSORIES CRYOMILL

02.480.0003



Autofill 150 l (recommended), incl. connection tube and safety valve

02.480.0002



Autofill 50 l, incl. connection tube and safety valve

05.871.0001



Connection tube, incl. safety valve (for LN2 supply provided by customer)

02.706.0304



Adapter for use of 2/4 grinding jars, 5 ml

02.706.0303



Adapter for use of 2/4/6 reaction vials, 2 ml

02.706.0360

Adapter for use of 2/4/6 reaction vials, 2 ml stainless steel tubes

22.749.0001



Safe-lock reaction vials 2 ml, 1000 pcs.

22.749.0008



Reaction vials made of stainless steel 316L, 2.0 ml, 10 pcs.

99.200.0016



IQ/OQ Documentation for CryoMill

22.111.0001

Gasket for grinding jar 5 ml, 10 pcs.

22.085.0019

Gasket for grinding jar 10 ml, 10 pcs.

22.085.0022

Gasket for grinding jar 25 ml, hardened steel or stainless steel, 10 pcs.

22.085.0023

Gasket for grinding jar 25 ml, zirconium oxide, 10 pcs.

22.085.0024

Gasket for grinding jar 35 ml, 10 pcs.

22.085.0025

Gasket for grinding jar 50 ml, 10 pcs.

GRINDING BALLS

HARDENED STEEL

05.368.0029  5 mm Ø

05.368.0030  7 mm Ø

05.368.0059  10 mm Ø

05.368.0032  12 mm Ø

05.368.0108  15 mm Ø

STAINLESS STEEL

05.368.0034  5 mm Ø

05.368.0035  7 mm Ø

05.368.0063  10 mm Ø

05.368.0037  12 mm Ø

05.368.0109  15 mm Ø

05.368.0062  20 mm Ø

05.368.0105  25 mm Ø

ZIRCONIUM OXIDE

05.368.0146  7 mm Ø

05.368.0094		10 mm Ø
05.368.0096		12 mm Ø
05.368.0113		15 mm Ø

PTFE WITH STEEL CORE

05.368.0045		10 mm Ø
05.368.0046		12 mm Ø
05.368.0114		15 mm Ø
05.368.0047		20 mm Ø