



颚式粉碎机 BB 400

符合人体工程学的操作和大进料尺寸


颚式破碎机**BB 400**是专为接受高达**220 x 90 mm**的大进料尺寸而设计的。间隙宽度可在**0-30 mm**范围内无级设置。根据样品材料的不同，可获得小于**2 mm**的颗粒尺寸，从而产生**50倍**的破碎比。**BB 400**的处理量可达**400 kg/h**。

RETSCH 颚式破碎机系列有 8 种不同规格。产量和最终细度取决于破碎机类型、所选间隙宽度和样品材料的破碎特性。

设计坚固、操作和清洁简单是 RETSCH 颚式破碎机型号的特点。少量样品可分批处理，大量样品可连续操作；

稳固 & 多功能

- | 前门便于清洁破碎室
- | 无限/连续间隙宽度设置
- | 适用于无污染研磨的多种材料
- | 带出口的收集装置，可用于连续破碎
- | 吸尘接头
- | 刻度显示鄂板磨损
- | 为特定应用选择料斗

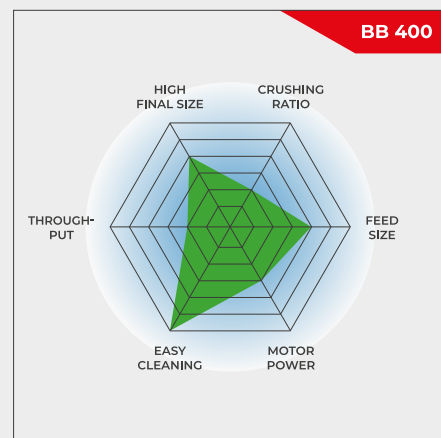


符合人体工程学的操作和大进料尺寸

颚式粉碎机 BB 400

优于细节

??10????????????????BB 400??BB400??
??BB 400??5????????????????????
????????BB400??

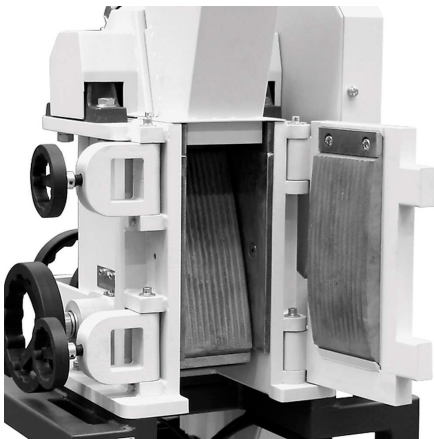

前置门便于清洁




除尘连接

BB 400鄂式破碎机配件

破碎鄂板

破碎鄂板和耐磨板有五种不同的材质可供选择，具体取决于样品和要进行的分析：

- | 锰钢
- | 不锈钢
- | 碳化钨
- | 防重金属研磨钢
- | NiHard4 钢铁

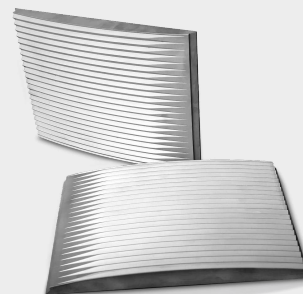

料斗类型

BB 400??
????????10??BB
400? ???



收集容器

????????10??
????????????26??



颚式粉碎机 BB 400

典型樣品材料

RETSCH强大的颚式破碎机非常适合对建筑材料、矿石、花岗岩、氧化物陶瓷、石英、矿渣、硅、煤、钨合金、水泥熟料等进行初步破碎。



水泥熔渣



陶瓷



催化剂




岩芯

为了找到解决方案，请访问我们的应用数据库。

颚式粉碎机 BB 400

作用原理

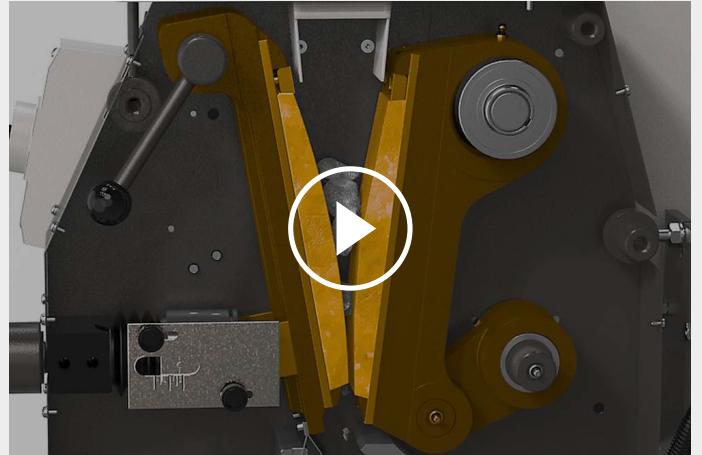[illegible]

[点击观看视频](#)

颚式粉碎仪 BB 400

作用原理

应用	粗粉碎
应用领域	化学 / 合成材料, 地质 / 冶金, 工程/电子, 建筑原料, 环境 / 资源回收利用, 玻璃/ 陶瓷
样品特征	中硬性, 硬的, 脆性的, 硬的
粉碎原理:	压力
最大进样尺寸	< 220 x 90 mm
最终出料粒度 *	< 2 mm
处理量*:	400 kg/h
研磨套件材料:	锰钢, 碳化钨, 防重金属研磨钢, NiHard4
颚板宽度	220 x 90 mm
研磨间隙宽度调节	0 - 30 mm
间隙宽度显示	规模
零点-校准	是
磨损率显示	是的, 可达负10毫米
可移动进料漏斗	是
加配吸尘器	是
中央润滑 (选项)	无
安装装备型供货	可选项
接收槽容积	10 l
底座	变形框架 (可调整为连续破碎)
驱动	三相电机
驱动功率	5.5 kW
电源数据:	不同电压
电源接头:	三相
防护类型	IP 55
机体尺寸 (宽x高x纵深)	695 x 1399 x 719 mm
净重	~ 400 kg
标准	CE

受样品材料性质和仪器配置/设定的影响

www.retsch.cn/bb400

ORDER DATA

JAW CRUSHER BB 400

Jaw Crusher BB 400 complete with breaking jaws, wearing plates, door, feed hopper, 10 l collecting receptacle and base frame

				粉碎颚板 / 刮磨板
21.002.0001		BB 400	3~ 400 V, 50 Hz	manganese steel / hardened steel
21.002.0012		BB 400	3~ 400 V, 50 Hz	NiHard 4 / hardened steel
21.002.0002		BB 400	3~ 400 V, 50 Hz	tungsten carbide / hardened steel
21.002.0004		BB 400	3~ 400 V, 50 Hz	heavy metal free steel / heavy metal free steel

其余可选的电源配置机型价格相同

CYCLONE FOR BB 400

CYCLONE SET WITH HOLDER, CONNECTING ELEMENT, LID WITH CONNECTOR FOR VACUUM CLEANER AND SAMPLE OUTLET
(PLEASE ORDER CASSETTE WITH OUTLET SEPERATELY)

22.935.0039	Cyclone set with collecting receptacle 30 l
03.183.0142	Cassette with outlet for use with cyclone
02.605.0420	Adjustable steel connection tubing (316L) for cyclone set
05.707.0208	FDA certified hose, 1.5 m
22.748.0017	Industrial vacuum cleaner Attix, 230 V, 50/60 Hz

ACCESSORIES JAW CRUSHER BB 400

BREAKING JAWS BB 400

23.048.0029	颚板, 锰钢, 1对
23.048.0045	颚板, NiHard 4, 1对
23.048.0030	颚板, 碳化钨, 1对
23.048.0031	颚板, heavy metal free steel, 1对

other materials available on request

WEARING PLATES BB 400

23.711.0021	摩擦板, 淬火钢, 1对
23.711.0020	摩擦板, 碳化钨, 1对
23.711.0019	摩擦板, heavy metal free steel, 1对

other materials available on request

FURTHER ACCESSORIES BB 400

23.003.0016	收集容器, standard steel, powder coated, 5 litres
02.054.0058	收集容器, standard steel, powder coated, 10 litres
23.003.0024	收集容器, 不锈钢, 26 litres complete with connection flange and carrier with wheels
02.785.0226	Collecting funnel with outlet for continuous use, powder coated
23.785.0018	Hopper for drill cores
23.003.0019	Infeed hopper, 不锈钢, 10 litres incl. connection rod and manual pusher
03.025.0137	Adapter for vacuum cleaner
01.657.0031	Raised access step