



颚式粉碎机 BB 500

功率和最终细度的最佳组合


颚式破碎机**BB 500**是专为接受高达**110毫米**的大进料尺寸而设计的。间隙宽度可在**0-11 mm**范围内无级设置。根据样品材料的不同，在一次工作中可获得小于**0.5 mm**的粒径，从而获得**220倍**的良好破碎比。**BB 500**的处理量可达**500 kg/h**。

RETSCH 颚式破碎机系列有 8 种不同规格。产量和最终细度取决于破碎机类型、所选间隙宽度和样品材料的破碎特性。

设计坚固、操作和清洁简单是 RETSCH 颚式破碎机型号的特点。少量样品可分批处理，大量样品可连续操作；

稳固 & 多功能

- | 极好的破碎率
- | 连续鄂板间隙设置
- | 适用于无污染研磨的多种材料
- | 吸尘接头
- | 带防溅保护的可拆卸料斗
- | 鄂板旋转180°可延长使用寿命

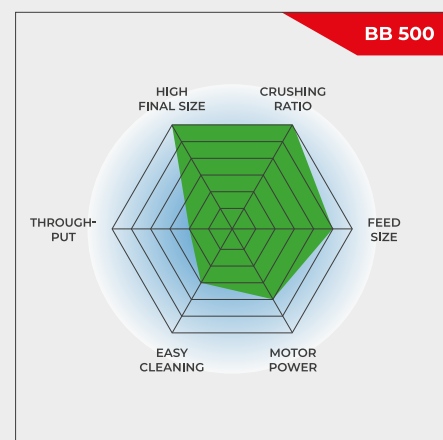

功率和最终细度的最佳组合


颚式粉碎机 BB 500

优于细节

BB 500 110 mm 0.5 mm 90% 15 l BB 500



用于连续间隙宽度设置的刻度

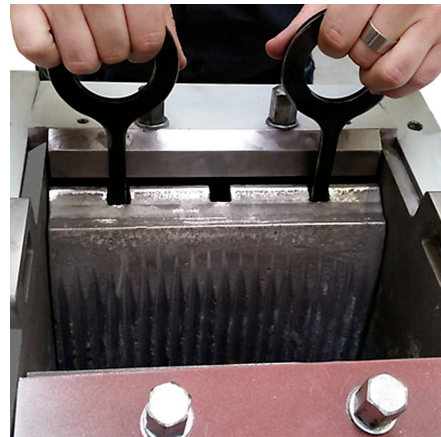

鄂板和耐磨板可以更换

生产线版本

BB 500



颚式粉碎机 BB 500

典型樣品材料

RETSCH强大的颚式破碎机非常适合对建筑材料、矿石、花岗岩、氧化物陶瓷、石英、矿渣、硅、煤、钨合金、水泥熟料等进行初步破碎。



熔铁



硅



石英岩




硅酸钙

为了找到解决方案，请访问我们的应用数据库。

颚式粉碎机 BB 500

作用原理

BB 500??
??该运动确保在冲程循环中保持一致
的间隙宽度，以便在一个工作步骤中将样品粉碎至设定的细
度。两个巨大的飞轮将高冲击力传递给破碎鄂板。由于鄂板
可以180°旋转，创新设计允许使用鄂板的两侧。????????????
??
????????????????

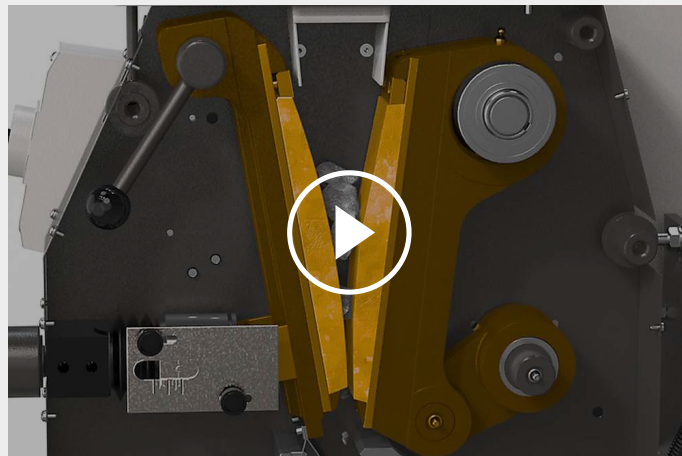

[点击观看视频](#)

颚式粉碎机 BB 500

技术参数

应用	粗粉碎
应用领域	化学 / 合成材料, 地质 / 冶金, 工程/电子, 建筑原料, 环境 / 资源回收利用, 玻璃/ 陶瓷
样品特征	中硬性, 硬的, 脆性的, 硬的
粉碎原理:	压力
最大进样尺寸	< 110 mm
最终出料粒度 *	90% < 0.5 mm
处理量*:	500 kg/h
研磨套件材料:	锰钢, NiHard 4
颚板宽度	250 x 130 mm
研磨间隙宽度调节	0 - 11 mm
间隙宽度显示	规模
可翻置漏斗	是
加配吸尘器	是
中央润滑 (选项)	无
安装装备型供货	可选项
接收槽容积	15 l
驱动	三相电机
驱动功率	7.5 kW
电源数据:	不同电压
电源接头:	三相
防护类型	IP 55
机体尺寸 (宽x高x纵深)	930 x 1400 x 1080 mm
净重	~ 1000 kg
标准	CE

受样品材料性质和仪器配置/设定的影响

www.retsch.cn/bb500

ORDER DATA

JAW CRUSHER BB 500

Jaw Crusher BB 500 complete with breaking jaws, wearing plates, feed hopper, collecting receptacle and base frame

				粉碎颚板 / 刮磨板
21.003.0001		BB 500	3~ 400 V, 50 Hz	manganese steel / hardened steel
21.003.0006		BB 500	3~ 400 V, 50 Hz	NiHard 4 / hardened steel

其余可选的电源配置机型价格相同

VERSION FOR INTEGRATION

Version for integration into a system at the customer's location, with manufacturer's declaration according to EC machine directive, with three-phase motor, without hopper, frame and protective motor switch

21.003.0003		BB 500	3~ 400 V, 50 Hz	manganese steel / hardened steel
21.003.0007		BB 500	3~ 400 V, 50 Hz	NiHard 4 / hardened steel

CYCLONE FOR BB 500

CYCLONE SET WITH HOLDER, FDA CERTIFIED HOSE, LID WITH CONNECTOR FOR VACUUM CLEANER AND SAMPLE OUTLET
(PLEASE ORDER CASSETTE WITH OUTLET SEPARATELY)

22.935.0039	Cyclone set with collecting receptacle 30 litres
03.183.0154	Cassette with outlet for use with cyclone
02.605.0420	Adjustable steel connection tubing (316L) for cyclone set
05.707.0208	FDA certified hose, 1.5 m
22.748.0017	Industrial vacuum cleaner Attix, 230 V, 50/60 Hz

ACCESSORIES JAW CRUSHER BB 500

BREAKING JAWS BB 500

23.048.0033



颚板, 锰钢, 1对

23.048.0036

颚板, NiHard 4, 1对

other materials available on request

WEARING PLATES BB 500

23.711.0025

摩擦板, 淬火钢, 1对

23.711.0030

摩擦板, NiHard 4, 1对

other materials available on request

FURTHER ACCESSORIES BB 500

23.003.0025

收集容器, standard steel, powder coated, 15 litres

05.045.0528

Collecting receptacle for cyclone set 22.935.0039, 30 litres